



Moving Health Forward

# Dyna-Form<sup>®</sup> Mercury Advance Trolley

User Manual





# Dyna-Form® Mercury Advance Trolley

The **Dyna-Form® Mercury Advance Trolley** Mattress has been designed especially for users considered to be that of **Very High Risk** of pressure ulcer development in potentially lengthy A&E waiting environments. This product can be made to measure in order to fit most trolley designs and is ideal for use in busy A&E departments that have an extensive patient turnover. Offering high levels of patient comfort, this system has the facility to “step up” to that of a dynamic mattress when clinically required. Similarly, the mattress function can be downgraded as the patient’s condition improves. These features make it particularly beneficial for use in an acute environment and help reduce moving and handling, logistical and decontamination impact costs. A higher maximum weight capacity, up to 40 stone / 254kg, allows the product to meet the modern challenges of those heavier clients. All component parts are interchangeable and replaceable, maximising product life and reducing environmental impact. The outer cover comprising a high frequency welded, multi stretch and vapour permeable fabric satisfies the strictest infection control policies.

## Contents

1. Introduction .....	4
2. Quick Reference Guide & Frequently Used Functions.....	4
3. Troubleshooting .....	6
4. Installation .....	7
5. Operation.....	7
6. Transportation.....	8
7. Warnings .....	8
8. Maintenance Procedures.....	9
9. Technical Specification .....	10
10. Technical Data.....	14
11. Optimum Conditions for Use .....	14
12. Symbols Guide & Contraindications for Use .....	14
13. Detachable / Removable Parts .....	15
14. Disposal .....	15

**Intelligent®**  
Pressure Care  
Management | Making a  
Measurable Difference

## 1. Introduction

The Mattress consists of a foam head cell and series of 14 transverse air cells, each containing a unique foam profiled insert, which are in turn held within a foam U Core, all protected by a vapour permeable waterproof cover. The single head end cell and the formers consist of foam only. The transverse cells are arranged into alternate pairs of A and B cells which are filled and emptied in sequence.

In Static Mode, the mattress attains the pressure reducing properties of the Dyna-Form Mercury static foam mattress (details available on request), whilst in Alternating Mode the mattress is able to offer similar properties to a pressure relieving dynamic system.


The digitally controlled Power Unit controls a pump that allows air to flow into, or out of the air cells as required according to the operating mode selected. It also maintains the air pressure within the mattress at the required level and controls the action of the audible/visual Audible Warning System in the event of mains supply failure or over or under inflation pressure. A CPR Valve located at the pump end of the umbilical hose permits the rapid deflation of the Mattress in an emergency.



## 2. Quick Reference Guide (Frequently used functions)

This is a quick reference guide for the **Dyna-Form Mercury Advance Trolley System**

**Product Code MAT1260001**

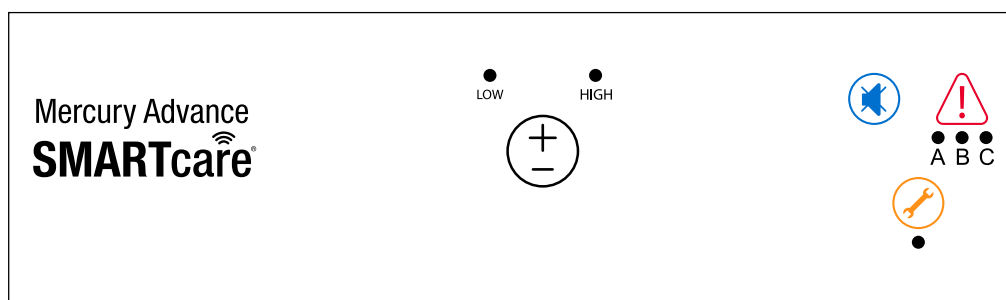

### Power Switch Audible Warning Reset

The power switch simply switches the mains power to the pump on and off.

When the pump detects an Audible Warning condition, this can be silenced (see page 5) and re-set by switching the pump off and then back on again.



### CPR Valve

Please ensure that the CPR connector is always placed fully home, prior to inflating the mattress. NB: The mattress will NOT inflate properly should this not be the case.

The CPR connector is only to be used in the event of a clinical emergency for priority use. However, disconnecting this function will cleverly deflate air rapidly from the mattress in readiness for transport / static mode.



## LED Mode Settings

This symbol when illuminated (The green indicator light) is not used to indicate that the equipment is on or ready for use.

When a patient requires a true dynamic function or indeed more pressure in the cells, as they may be uncomfortable or feel as though the support surface is too soft or unstable, then please select a “High” setting (pressure 26mmHg). This must only be used by a trained clinician as often too high pressures can further agitate certain patient conditions.

When a patient requires less pressure in the cells, as they may be uncomfortable or indeed hyper sensitive to cell movement or indeed if the patient is still reddening further, then please select a “Low” setting (pressure 18mmHg). This must only be used by a trained clinician.

This function is used to silence the Audible Warning. The LED will remain lit if the Audible Warning has been silenced previously, however a fault is still detected. Refer to the power switch in order to re-set fully. If the Audible Warning continues to sound repeatedly, along with an illuminated light, then an engineer must be called.

This symbol indicates an “Audible Warning Failure”. Please see troubleshooting guide below for how to re-set.

This symbol when illuminated indicates a Service is required. DHG recommends an annual service. The service light will illuminate every 8760 hours of operation (one continuous year running).

Note: Please ensure (when available) that all securing straps on the base of the mattress are secured onto the MOVING PARTS of the bed frame.

For shut down procedure, see 4.2 Power Unit (Pump) section.

● LOW      ● HIGH



Power On / Off  
True Dynamic / Firmer Setting

● LOW      ● HIGH



Low / Comfort Pressure Setting



Silence Audible Warning









● ● ●  
A B C

Audible Warning Failure



Service Indicator

### 3. Troubleshooting

Warning / Fault	Cause	Solution
<b>Control Unit does not operate; no display lights illuminate</b>	The Control Unit may not be attached to a power source or a fuse may need replacing.	<ol style="list-style-type: none"> <li>1. Check the Control Unit is connected to mains power outlet with the correct voltage.</li> <li>2. Check the Control Unit is switched on. Switch off and unplug the unit before restarting.</li> <li>3. Check the mains plug fuse (5 AMP) then check both Control Unit fuses (1 AMP) – fuses can be released using a screwdriver to push and turn.</li> </ol> <p>⚠ Do not try to open the Control Unit. Opening the unit could cause personal injury or equipment damage.</p> <p>⚠ Ensure the replacement of fuses is carried out accordance with local legislation.</p>
<b>Warning LED C + audible warning</b> 	Mains failure / Other See above, plus:	<ol style="list-style-type: none"> <li>1. Reset the audible warning – turn off power and press the audible warning mute button.</li> <li>2. Check the handle is intact, ensuring all four sealing connectors are firmly fitted to the control unit and the air hoses. Check the CPR tag is attached and all sealing connectors are firmly secure.</li> <li>3. Check all air hoses along the inside of the mattress – each should be firmly connected. Check each air cell is securely attached to its connecting air pipe.</li> <li>4. Check all cells, pipes and hoses for any air leakage.</li> <li>5. Switch on power.</li> </ol>
<b>Warning LED B + audible warning</b> 	Pressure too low	<ol style="list-style-type: none"> <li>1. Reset the warning – turn off power and press the audible warning mute button.</li> <li>2. Check the handle is intact, ensuring all four sealing connectors are firmly fitted to the control unit and the air hoses. Check the CPR tag is attached and all sealing connectors are firmly secure.</li> <li>3. Check all air hoses along the inside of the mattress – each should be firmly connected. Check each air cell is securely attached to its connecting air pipe.</li> <li>4. Check all cells, pipes and hoses for any air leakage.</li> <li>5. Check that the air filter cover is correctly secured and the air filter is clean.</li> <li>6. Switch on power.</li> </ol>
<b>Warning LED B+C + audible warning</b> 	Pressure too low / Air pipe kinked	<ol style="list-style-type: none"> <li>1. Check all air hoses along the inside of the mattress – each should be firmly connected. Check each air cell is securely attached to its connecting air pipe.</li> </ol>
<b>Warning LED A + audible warning</b> 	Pressure too high	<ol style="list-style-type: none"> <li>1. Reset the warning – turn off power and press the audible warning button.</li> <li>2. Disconnect the air hoses to reduce pressure, reconnect when pressure has decreased.</li> <li>3. Check for twists in the air hoses between Mattress and Control Unit.</li> <li>4. Switch on power.</li> </ol>
<b>Warning LED A+B + audible warning</b> 	Alternating Mode Failure (no alternation)	<ol style="list-style-type: none"> <li>1. Reset the warning – turn off Power and press the audible warning mute button.</li> <li>2. Disconnect the air hoses to reduce pressure – reconnect when pressure has decreased.</li> </ol>
<b>Warning LED A,B+C + audible warning</b> 	Initialising Failure	<ol style="list-style-type: none"> <li>1. Press the audible warning mute button to silence the audible warning.</li> <li>2. Check the power cable is firmly plugged into the mains outlet and the Control Unit; and check the mains power is switched on.</li> <li>3. Check the Control Unit fuse (1 AMP) – fuses can be released using a screwdriver to push and turn.</li> </ol>

## 4. Installation

### 4.1 Mattress (This is the applied part type BF)

Place the Dyna-Form Mercury Advance Trolley Mattress directly on to the bed platform ensuring that the Blue multi-stretch waterproof cover is on top and that the umbilical hose is located at the left hand corner at the foot end of the bed. Note: The umbilical hose can be located inside the cover under the “Open Here for Air Inlet” printed in the bottom left hand corner of the mattress.

Wipe the mattress down before covering the mattress with a loose-fitting sheet.

#### Static Mattress Use

The Dyna-Form Mercury Advance Trolley Mattress can be used as a pressure reducing mattress for patients at High Risk of pressure ulcer damage without the need to attach the pump.

#### Alternating Mattress Use

If / When required, the Dyna-Form Mercury Advance Trolley Mattress can be used as an alternating mattress system by attaching the Dyna-Form Mercury Advance Trolley Power Unit (Pump) system.

**No other system should be attached to the mattress as the design settings and internal air pressure properties of the Dyna-Form Mercury Advance Trolley Power Unit (Pump) are specific to this mattress only.**

The Dyna-Form Mercury Advance Trolley is a replacement mattress system and should NOT be placed on top of any existing mattress.

The startup time from static to dynamic mode is immediate.

### 4.2 Power Unit (Pump)

Hang the Power Unit (Pump) onto the footboard. The mounting hooks swivel to suit the thickness of the footboard or rail. Connecting the Umbilical Hose to the Power Unit (Pump), place the 3-pin electrical plug into the wall outlet and switch on:

- (a) Open the zip located at the bottom left hand side of the mattress and pull out the Blue Umbilical hose.
- (b) Attach the Blue Umbilical Hose to the Power Unit (Pump) by connecting the air connector at the end of the Umbilical Hose to the air inlet connector at the bottom left hand side of the pump. Ensure that the Red CPR Release button is located on top of the Air Inlet connector after connection is complete.
- (c) Re-close the zip as far as possible without clamping the Blue Umbilical Hose to ensure the mattress and air cells are sealed within the cover.
- (d) Shut down is the reverse of items a, b & c above.

## 5. Operation

Attach the mains cable to the pump by inserting the “kettle” type connector into the recess located on the left hand side of the pump. The mains cable has been designed specifically as a removable part to aid in easy replacement should it become damaged in use.

The mains plug should be turned off and removed from wall socket as a means of isolation.

Plug the mains cable into a suitable 230v mains socket and switch on the Power Unit using the on/off switch.


After the pump has been turned on both the “High” and the “Low” lights will flash together intermittently until the pump has attained its initial operating pressure. Once the pump has attained its initial operating pressure the “Low” light will stay on constantly and the mattress is ready for use.



### 5.1 Low / High Settings

The Dyna-Form Mercury Advance Trolley Mattress, in Alternating Mode, has two pressure settings. The initial setting that the pump will revert to upon set up is “Low”. The “Low” comfort setting is ideal for the lighter patient or those who feel discomfort when on a normal alternating air type mattresses system. However, for patients with existing pressure damage or those at Very High Risk, it is recommended that dependant on the clinical judgement of the clinician, the “High” setting is activated by pressing the +/- button once, which is located on top of the pump.

In “High” Mode the pump attains more of the characteristics of an alternating air mattress system whilst still utilising the advantages of the static foam inserts. Repeatedly pressing the ‘mode’ button enables the Low & High modes to be selected in turn.

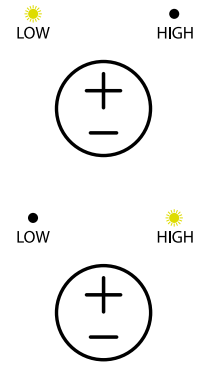

### 5.2 CPR Deflation

The CPR system consists of a manually operated button located on the Air Inlet connector attached to the pump. By pressing the Red Button, which will release the connector locking system, the user can remove the connector unit which will deflate the mattress air system back to that of a static foam mattress.

Note: After a short period as the Mattress deflates the ‘Low Pressure’ Audible Warning is activated and can be cancelled by switching the Power Unit off.



### 5.3 Troubleshooting

For assistance (if needed) in setting up, using or maintaining the Mercury Advance System, or to report unexpected operation or events, please contact Direct Healthcare Group on the contact details on the reverse of this manual.

## 6. Transportation

To change the location of the mattress, remove the Umbilical Cord and allow the mattress to return to its Static Mattress form. Switch off the Power Unit (Pump) using the on/off switch and disconnect the electrical supply cable from the mains socket. The mattress can now be moved to a new location where it must immediately be reconnected to the mains electrical supply and the Power Unit (Pump) switched back on. Once the Mattress has been refilled, the ‘Alternating’ mode will automatically revert back to the Low setting and should be reselected to High should this be desired by the clinician.

Warning: The Mattress will not ‘alternate’ when disconnected from the Power Unit (Pump) and/or the mains electrical. Also refer to environmental conditions section at rear of this manual.

## 7. Warnings

Warning conditions are indicated by a flashing red display accompanied by an audible warning.

In each case the user should respond by turning the Power Unit’s switch off and investigating the cause.


### 7.1 High Pressure Warning

This condition could be caused, for example by a kinked Umbilical Hose or visitors, and others, sitting suddenly on the Mattress.



### 7.2 Low Pressure Warning

This condition could be caused, for example, by incorrect fitting of the air inlet connector, opening of the CPR Valve or a leak in the Mattress due to a cut or puncture.

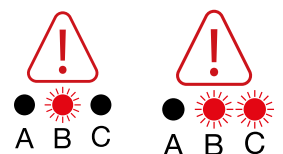

### 7.3 Mains Failure Warning

This condition may be caused, for example if Mains power is lost.





### 7.4 Alternating Mode Failure (no alternation)

This will be indicated by a warning LED on A and B and an audible warning.

1. Reset the warning – turn off Power and press the audible warning mute button.
2. Disconnect the air hoses to reduce pressure – reconnect when pressure has decreased.



### 7.5 Initialising Failure

This will be indicated by a warning LED on A, B and C and an audible warning.

1. Press the audible warning mute button to silence the audible warning.
2. Check the power cable is firmly plugged into the mains outlet and the Control Unit; and check the mains power is switched on.
3. Check the Control Unit fuse (1 AMP) – fuses can be released using a screwdriver to push and turn.



## 8. Maintenance Procedures

### 8.1 Safety Warning

Only qualified technicians trained or formally approved by Direct Healthcare Group Ltd. in the operation and maintenance of Direct Healthcare Group products may carry out maintenance, modification or repair work on the equipment. Unqualified personnel attempting to work on Direct Healthcare Group Power Units risk serious injury to themselves and others and possibly death by electrocution. Inlet fuse NOT to be replaced by operator or patient, to be replaced by service personnel only.

**Warning – Do not modify this equipment without authorisation of Direct Healthcare Group.**

#### 8.1.1 Servicing

Direct Healthcare Group recommend that the control unit be serviced annually from installation. The service light will illuminate after 8760 hours of operation (one year of continuous operation). The unit contains no user serviceable parts and should only be carried out by persons as described in section 8.1. DHG will make available on request service manuals, component parts lists and other information necessary for any suitably qualified person (As in 8.1) to carry out repair or service the system. For Service, maintenance and any questions regarding this please contact DHG.

### 8.2 Cleaning Procedures

**Warning: Before cleaning the System make sure that the Power Unit (Pump) is disconnected from the mains electricity supply.**

**Do not** immerse the Power Unit (Pump) in water or other fluids.

**Do not** autoclave, nor use phenol for cleaning.

Do wash hands before commencing the cleaning process.

Wear appropriate protective clothing such as gloves, apron and a mask.

Ensure all work surfaces are cleaned before and after contact with the Mattress.

### 8.3 Warning – Cleaning the Mattress

1. Cleaning should take place before and after use, and between patients.
2. With cover left on the Mattress disconnect the Mattress from the Power Unit (Pump).
3. Clean the surface of the wash down table with Hypochlorite solution or equivalent disinfectant.
4. Wash Mattress top using hot water (60 degrees C) containing detergent – dry with a paper towel.
5. Use a Hypochlorite solution 1,000 parts per million available chlorine. For heavy contamination use a Hypochlorite solution 10,000 parts per million available chlorine. Please ensure thorough rinsing after cleaning.
6. Using suitable brush, hot water, detergent or Hypochlorite solution, clean Umbilical Hose and CPR Valve. Dry with paper towel.
7. If required, the Mattress Cover may be removed and machine-washed at a temperature of 80 degrees C, for not less than 10 minutes. The individual Air Cells can be wiped down with established disinfectants.
8. To avoid shrinkage of the cover line dry in an indoor clean environment or tumble dry on a low heat setting not exceeding 40 degrees C and not for longer than 10 minutes. Covers must be thoroughly dried before re-fitting to the mattress.

### 8.4 Warning – Cleaning the Power Unit (Pump)

The Power Unit can be cleaned by wiping with a cloth dampened with a detergent solution or Hypochlorite solution. Also refer to symbol chart.

#### 8.4.1 Warning

Ensure the Mercury Advance System is not exposed to:

1. Excessive heat sources e.g. fires, radiators etc.
2. Water, particularly immersion of the pump.

## 9. Technical Specification

### Declaration – Electromagnetic Emissions - for all ME EQUIPMENT and ME SYSTEMS

<p>Guidance and manufacturer’s declaration – Electromagnetic Emission</p> <p>The MAT1260001 is intended for use in the electromagnetic environment specified below.</p> <p>The customer or the user of the system should ensure that it is used in such an environment.</p>		
Emission Test	Compliance	Electromagnetic Environment – Guidance
RF emissions CISPR 11	Group 1	The system uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emission CISPR 11	Class B	The system is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

## 9. Technical Specification

### Declaration – Electromagnetic Immunity

Guidance and manufacturer’s declaration – Electromagnetic Immunity

The MAT1260001 is intended for use in the electromagnetic environment specified below.

The customer or the user of the system should ensure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level		Compliance Level		Electromagnetic Environment – Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air		±6 kV contact ±8 kV air		Floors should be wood, concrete or ceramic tile. If floor are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output line(s)		±2 kV for power supply lines ±1 kV for input/output line(s)		Mains power quality should be that of a typical commercial or hospital environment.
Surge Immunity Test IEC 61000-4-5	± 1 kV line(s) to line(s)		±1 kV differential mode		Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Voltage Dips %U <sub>T</sub>	Period (Cycles)	Voltage Dips %U <sub>T</sub>	Period (Cycles)	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Span system requires continued operation during power mains interruptions, it is recommended that the system be powered from an uninterruptible power supply or a battery.
	30	25	30	25	
	60	5	60	5	
	>95	0.5	>95	0.5	
	Voltage Interruption % U <sub>T</sub>	Seconds	Voltage Interruption % U <sub>T</sub>	Seconds	
	>95	5	>95	5	
Power frequency (50Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m		Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.	

NOTE U<sub>T</sub> is the a.c. mains voltage prior to application of the test level.


## 9. Technical Specification

### Declaration – Electromagnetic Immunity – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacturer’s declaration – Electromagnetic Immunity

The MAT1260001 is intended for use in the electromagnetic environment specified below.

The customer or the user of the system should ensure that it is used in such an environment.

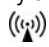Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF IEC 61000- 4-6	3 V <sub>rms</sub> 150 kHz to 80 MHz	3 V <sub>rms</sub>	Portable and mobile RF communications equipment should be used no closer to any part of the CT515, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 61000-4-3	3 V/m (Professional Healthcare Environment) 10 V/m (Home Healthcare Environment) 80 MHz at 2.7 GHz	10 V/m	Recommended separation distance $d = 1.167\sqrt{P}$ $d = 1.167\sqrt{P}$ 80 MHz to 800 MHz $d = 2.333\sqrt{P}$ 800 MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, <sup>a</sup> should be less than the compliance level in each frequency range. <sup>b</sup> Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Span system is used exceeds the applicable RF compliance level above, the system should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the system.

b. Over the frequency range 150kHz to 80MHz, field strengths should be less than 3V/m.

## 9. Technical Specification

### Recommended separation distances between portable and mobile RF communications equipment and the EQUIPMENT or SYSTEM – for ME EQUIPMENT or ME SYSTEM that are not LIFE – SUPPORTING

Recommended separation distances between portable and mobile RF communications equipment and the MAT1260001 Alternating Control Unit.			
The MAT1260001 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the system can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the system as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 KHz to 80 MHz $d = 1.167\sqrt{P}$	80 MHz to 800 MHz $d = 1.167\sqrt{P}$	800 MHz to 2.5 GHz $d = 2.333\sqrt{P}$
0.01	0.117	0.117	0.233
0.1	0.369	0.369	0.738
1	1.167	1.167	2.333
10	3.689	3.689	7.379
100	11.667	11.667	23.333
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

## 10. Technical Data

### 10.1 Power Unit (Pump)

Serial Number .....	As per label on rear of pump
Electrical Supply .....	220 – 240 volt, 50 Hz
Power Consumption .....	10 watts
Fuses .....	TA1H 250V
Protection against shock .....	Class 2
Noise Level .....	Approx. 30 dB (A)
Dimensions .....	245 x 160 x 95 mm
Weight.....	1.7 kg
Service Interval .....	12 months / 8760 hours
Expected life .....	5 years
Shelf life of parts .....	5 years

### 10.2 Mattress

Serial Number .....	Label on inside of mattress cover
Number of Air Cells .....	14 Air Cells / 1 Static Foam Cell
Dimensions .....	Made to Measure
Weight.....	10kg
Expected life of Mattress .....	5 years
Shelf life of Mattress parts .....	5 years

## 11. Optimum Conditions

(Applies to Mattress and Pump)

### 11.1 Environment Conditions for Use










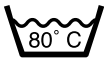


Transport.....	-25°C – +70°C
Storage .....	-25°C – +70°C
Usage .....	+5°C – +40°C
Humidity.....	10% – 93%
Atmospheric Pressure .....	700hPa – 1060hPa
Operational Altitude .....	≤ 2000m

### 11.2 Exposure

Exposure to direct sunlight, dust, lint and general debris is not considered to be an issue with the Mercury Advance System.

## 12. Symbols Guide

### Mattress Symbols

 WASH AT 80°	 TUMBLE DRY ON LOW	 DO NOT DRY CLEAN
 REFER TO USER MANUAL	 1434	 TYPE BF APPLIED PART
 DO NOT BLEACH	 DO NOT IRON	 NO SMOKING
 MAXIMUM USER WEIGHT LIMIT 254 KG / 40 STONES	 DO NOT USE SHARP INSTRUMENTS	 DO NOT USE PHENOL






#### WARNING

THIS IS A STATEMENT THAT ALERTS THE USER TO THE POSSIBILITY OF SERIOUS INJURY OR OTHERWISE ADVERSE REACTIONS WITH THE USE OR MISUSE OF THE DEVICE






#### CAUTION

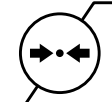THIS IS A STATEMENT THAT ALERTS THE USER TO THE POSSIBILITY OF A PROBLEM WITH THE SYSTEM ASSOCIATED WITH ITS USE OR MISUSE

### General Symbols

 CAUTION	 PROTECT FROM HEAT AND RADIOACTIVE SOURCES	 TEMPERATURE LIMITATION
 HUMIDITY LIMITATION	 ATMOSPHERIC PRESSURE LIMITATION	

### Pump (Unit) Symbols

 SERVICE REQUIRED	 REFER TO USER MANUAL	 DO NOT DISPOSE OF WITH HOUSEHOLD WASTE. PLEASE REFER TO DHS WEBSITE
 DOUBLE INSULATED CLASS II	 1434	

## IP21

IP: INGRESS PROTECTION  
2: PROTECTION AGAINST FINGERS OR OTHER OBJECT NOT GREATER THAN 80MM IN LENGTH AND 12MM IN DIAMETER  
1: PROTECTION FROM VERTICALLY DRIPPING WATER

**Contraindications For Use (Warning)**

The Mercury Advance System should not be used for patients with unstable fractures, gross oedema, burns, or intolerance to motion.

**General Information (Caution) (Warning)**

- There are no special skills required to operate the system.
- The Medical Professional is responsible for applying his/her best medical judgment when using the system.
- The electricity supply is of the type indicated on the Power Unit (pump).
- Check the mains lead is free from damage and is positioned so as not to cause an obstruction, or injury. E.g. Strangulation of a child or trip hazard.
- Ensure the mains lead cannot become trapped or crushed, e.g. by raising or lowering of the bed or bed rails or any other moving object.
- The power unit (pump) must only be used with a suitably approved power cord and plug set as supplied by DHG.
- The system is not to be used in the presence of flammable anaesthetics.
- Suitable for continuous use.
- Not suitable for sterilisation.
- Do not position the power unit to make it difficult to disconnect the power supply or plug.
- Do not place the System on or close to a source of heat ( i.e. a radiator) as this excessive exposure to heat may weaken the cell material.
- Do not use with hot water bottles or electric blankets.
- DHG strongly advise against smoking whilst the Power Unit (pump) is in use. This is to prevent accidental secondary ignition of items which may be flammable e.g. bed linen. The materials used in the manufacture of the Mercury Advance System comply with the required fire safety regulations.
- Do not use sharp objects on or near the mattress system as this will cause damage.
- Do not store in damp conditions.
- Do not use in an oxygen enriched environment.
- Not suitable for use in an Outdoor Environment.
- Intended for both Home Healthcare and Professional Healthcare environments.

- Do not connect to any other medical device or equipment.
- Correct fuse rating **MUST** be used. Failure to do so could result in the risk of a fire.
- The System should be cleaned after use or between patients. Refer to Cleaning section.
- All internal and external hoses must be free of twists, kinks. The external hose should also be properly connected and positioned so that the risk of obstruction or injury is eliminated.
- Do not use bleach, phenols. Chlorine based products which exceed 1000ppm. Solvents or alcohol based cleaners.
- All the above warnings and cautions together with safety considerations should be observed at ALL times during its use.
- Select correct setting 'Hi' or 'Low' as required. Care should be taken not to accidentally change settings once set. This may affect the desired requirement of the therapy. This could also be caused by pets, pests or children.
- This device does not emit radiation.

**13. Detachable/Removable Parts**

1. Mattress (Detached from the pump by removing the CPR connector). Part No. MAT1260002 (or variants of for the size)
2. Electric power cable. (Removed from the pump by pulling the cable away from the mains inlet on the side of the pump). Part No. SP021016

N.B. The battery is an integral part of the PCB and is not removable or changeable.

**Caution**

**Use of detachable parts not listed is not recommended by Direct Healthcare Group.**

**14. Disposal**

Please refer to DHG website for recommendations and responsibilities for disposal within the UK WEEE guidelines.



Moving Health Forward

**PRESSURE ULCER PREVENTION  
SAFE MOVING & HANDLING  
SPECIALIST THERAPIES  
BATHROOM SAFETY  
RENTAL & SERVICE SOLUTIONS**



**Direct Healthcare Group Ltd**  
Withey Court, Western Industrial Estate  
Caerphilly, United Kingdom  
CF83 1BF  
T: +44 (0) 800 043 0881  
F: +44 (0) 845 459 9832  
E: [info@directhealthcaregroup.com](mailto:info@directhealthcaregroup.com)



**EU Authorised Representative, Importer  
and Distributor:**

Direct Healthcare Group Sverige AB,  
Torshamnsgatan 35, SE-164 40, Kista  
Sweden  
T: +46 (0) 8 557 200

**Australian Sponsor:**

Direct Healthcare Group Pty Ltd.,  
68 Howe Street, Osborne Park,  
Western Australia 6017  
T: +61 (0) 423 852 810



LIT-00008 Issue 7  
Date: February 2022

[DIRECTHEALTHCAREGROUP.COM](http://DIRECTHEALTHCAREGROUP.COM)

