


Patient's name: \_\_\_\_\_

Social insurance number: \_\_\_\_\_

HygieneSling is easy to apply and facilitates dressing/undressing of the patient in connection with visits to the toilet. It features leg sections, which means that the patient does not sit on the sling. Instead, the sling provides support from the tailbone up over the head and under the legs. HygieneSling is designed for most lifting situations that involve dressing and undressing of the patient. HygieneSling is available in three models: without belt, with belt and a model with a belt and extra back support.

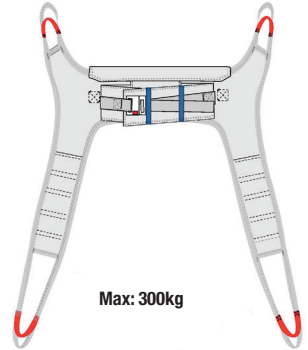

### Functional inspection

#### Visual and mechanical inspection

Check the condition and function of the sling regularly. Always inspect the product after laundering. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. Apply load to the device and check to ensure that clasps, handles, etc. withstand heavy load. If there are signs of wear, the product must be discarded.

#### Which size is used?

XX-Small     X-Small     Small     Medium     Large     X-Large     XX-Large

#### Which material is used?

Polyester

#### Accessories

SoftLegSupport     AdjustmentLoop     StrapsPadding

#### Which slingbar is used?

---

#### Leg support placement:



Crossed



Overlapped

## Instructions for use

Can the sling be left under the patient in the wheelchair?  Yes  No

Can the sling be used for transfers to/from the toilet?  Yes  No

Can the sling be used in showering situations?  Yes  No

Loop alternatives/AdjustmentLoop: \_\_\_\_\_

### Have staff that use the assistive device:

Read the instructions?  Yes  No

Received practical instruction in patient transfers?  Yes  No

### Instructions to personnel:

e.g.; Raise the bed side rail on the side towards which the patient is to be turned. Placement of head/arms/hands?  
Placement of legs/feet? Use bed functions.

---

---

---

---

Number of caregivers: \_\_\_\_\_

### Instructions to patient:

e.g.; In what way can the patient participate actively in the transfer situation? Which instructions are used?

---

---

---

---

### Other:

---

---

### Trial-fitted by:

Date: \_\_\_\_\_

Phone: \_\_\_\_\_

Phone: \_\_\_\_\_

### In the events of changes or questions, contact the trial-fitter.

Manuals and instructions for functional testing are available on DHG's website:

DHG-HEALTHCARE.COM

SystemRoMedic™ is a complete solution for easy transfers.  
The concept encompasses four product areas that meet  
most needs in different transfer situations.

**Transfer**, between two locations. **Positioning**, within the same location.  
**Support**, for mobility. **Lifting**, both manual and mechanical.

**DHG**  
Moving Health Forward

Withey Court, Western Industrial Estate  
Caerphilly, United Kingdom CF83 1BF  
T: +44 (0) 800 043 0881  
F: +44 (0) 845 459 9832  
E: info@directhealthcaregroup.com  
DHG-HEALTHCARE.COM