WendyLett2Way

SystemRoMedic™

WendyLett4Way


SystemRoMedic™

- Instructions for use – English
- Bruksanvisning – Svenska
- Brukermanual – Norsk
- Brugsvejledning – Dansk
- Käyttöohje – Suomi
- Gebrauchsanweisung – Deutsch
- Handleiding – Nederlands
- Manuel d'utilisation – Français
- Manuale utente – Italiano
- Manual de usuario – Español
- Návod k použití – Český
- Návod na použitie – Slovensky

WendyLett2Way

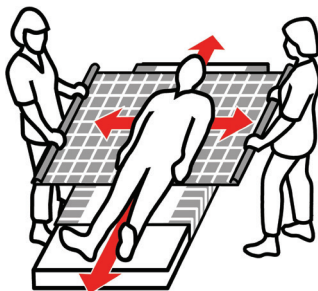WendyLett4Way

REF

1620
1631
1636
1636CH

REF

1637
1647
1648
1649

1652
1657
1637CH
1657CH

1687
1692

SWL: 300 kg/660 lbs
SWL 1649, 1657 & 1687: 200 kg/440 lbs

Table of contents

Product information	3
English	5
Svenska	8
Norsk	11
Dansk	14
Suomi	17
Deutsch	20
Nederlands	23
Français	26
Italiano	29
Español	32
Český	35
Slovensky	38
Symbols	41

**Product information/ Produktinformation/ Produktinformasjon/ Produktinformation/
Tuotetiedot/ Produktinformation/ Productinformatie/ Information produit/
Informazioni sul prodotto/ Información del Producto / Informace o výrobku / Informácie
o produkte**

Art. no.	Product/Product/ Tuote/ Produktname/ Product/ Produit/ Prodotto/ Producto / Výrobek	Size/ Storlek/ Størrelse/ Koko/ Größe/ Maat/ Taille/ /Dimensione/ Tamaño cm/ inch / Velikost	Colour/Färg/Farge/ Farve/ Väri/ Farbe/ Kleur/ Couleur/ Co- lore/ Color / Barva
----------	--	---	---

WendyLett - sliding sheet/ gliidlakan/ gliidelaken/ gliidelagen/ liukuvan telttakankaan/ piirtää arkki/ gleitlaken/ glijlaken/ feuille de glissement/ telone scorrevole/ sábana de deslizamiento / prostěradlo pro posouvání / posuvný list
(Can be combined with WendyLett2Way drawsheet or WendyLett4Way sliding drawsheet)

1609	WendyLett, 60 cm/24" sliding surface (bed 80/90 cm - 32"/36")	100 x 200/39 x 79	white
1611	WendyLett, 50 cm/20" sliding surface (bed 80/90 cm - 32"/36")	100 x 180/39 x 71	grey/white striped
1612	WendyLett, 50 cm/20" sliding surface (bed 80/90 cm - 32"/36")	140 x 180/55 x 71	grey/white striped
1613	WendyLett, 60 cm/24" sliding surface (bed 80/90 cm - 32"/36")	200 x 80/79 x 31	grey/white striped
1614	WendyLett, 60 cm/24" sliding surface (bed 80/90 cm - 32"/36")	140 x 200/55 x 79	white
1621	WendyLett, 75 cm/29" sliding surface (bed 80/90 cm - 32"/36")	140 x 205/55 x 81	grey/white striped
1622	WendyLett, 75 cm/29" sliding surface (bed 80/90 cm - 32"/36")	200 x 205/79 x 81	grey/white striped
1629	WendyLett, 60 cm/24" sliding surface (bed 80/90 cm - 32"/36")	100 x 200/39 x 79	grey/white striped
1634	WendyLett, 60 cm/24" sliding surface (bed 80/90 cm - 32"/36")	140 x 200/55 x 79	grey/white striped
1639	WendyLett, 70 cm/28" sliding surface (bed 105 cm/42")	200 x 200/79 x 79	grey/white striped
1667	WendyLett, 100 cm/40" sliding surface (bed 120 cm/48")	100 x 220/39 x 88	grey/white striped
1668	WendyLett, 100 cm/40" sliding surface (bed 120 cm/48")	140 x 220/55 x 88	grey/white striped
1669	WendyLett, 100 cm/40" sliding surface (bed 120 cm/48")	200 x 220/79 x 88	grey/white striped
1641	WendyLett, 70 cm/28" sliding surface (bed 90 cm/36"), tailored with elastic string	200 x 90/79 x 36	grey/white striped
1651	WendyLett, 80 cm/32" sliding surface (bed 105 cm/42"), tailored with elastic string	200 x 105/79 x 42	grey/white striped
1661	WendyLett, 95 cm/37" sliding surface (bed 120 cm/48"), tailored with elastic string	200 x 120/79 x 47	grey/white striped
1671	WendyLett, 110 cm/43" sliding surface (bed 140 cm/56"), tailored with elastic string	200 x 140/79 x 55	grey/white striped
1681	Wendylett, 70 cm/28" sliding surface	200 x 365/79 x 144	grey/white striped
1691	Wendylett, 70 cm/28" sliding surface	215 x 90/85 x 35	grey/white striped

WendyLett2Way - drawsheet/ draglakan/ stikkklakenet/ træklagen/ piirtää arkki/ zuglaken/ treklaken/ dessiner feuille/ foglio di copertura/ sábana travesera / napinací prostěradlo / žrebovací list

(To be combined with WendyLett sliding sheet)

1620	WendyLett2Way	90 x 180/35 x 70	white
1631	WendyLett2Way	100 x 200/39 x 79	white/grey striped
1636	WendyLett2Way	140 x 200/55 x 79	white/grey striped

WendyLett4Way - drawsheet/ draglakan/ stikkklakenet/ træklagen/ piirtää arkki/ zuglaken/ treklaken/ dessiner feuille/ foglio di copertura/ sábana travesera / napinací prostěradlo / žrebovací list

(To be combined with WendyLett sliding sheet)

1637	WendyLett4Way	140 x 200/55 x 79	grey/white checked
1647	WendyLett4Way	200 x 200/79 x 79	grey/white checked
1648	WendyLett4Way, inco	75 x 200/30 x 79	grey (inco) checked
1649	WendyLett4Way, inco with handles	140 x 120/55 x 47	grey (inco)
1657	WendyLett4Way, with handles	140 x 190/55 x 75	grey/white checked
1652	WendyLett4Way	140 x 260/55 x 102	grey/white checked
1687	WendyLett4Way, with handles	200 x 190/79 x 75	grey/white checked
1692	WendyLett4Way	215 x 200/85 x 79	grey/white checked

WendyLean - pillowcase with handles/ örngott med handtag/ putetreck med håndtak/ hovedpudebetræk med håndtag/ tynnyliina kahvat/ Kissenbezug mit Zugschlaufen/ kussenloop met handgrepen/ Taie d'oreiller avec des poignées/ federa con maniglie/ unda de almohada con asas / povlak na polštář s úchyty / obliečka na vankúš s rúčkami

1680	WendyLean	50 x 60/20 x 24	grey/white striped
------	-----------	-----------------	--------------------

Product combinations, Produktkombinationer, Produkt kombinasjon, Produkt kombination, Tuoteyhdistelmä, Produktkombination, Productcombinatie, Combinaison de produit, Combinazione di prodotti, Combinación del producto, Kombinace výrobků, Kombinácie produktov,

Art. no. Combination

1610	1609+1620
1615	1614+1620
1630	1629+1631
1635	1634+1636
1640	1639+1636
1642	1641+1636
1638	1634+1637
1628	1634+1657
1658	1634+1649
1643	1639+1637
1644	1641+1637
1645	1641+1647
1646	1641+1657
1650	1669+1636
1656	1651+1657
1655	1651+1637
1658	1634+1649
1665	1669+1647
1699	1639+1649

WendyLett2Way

SystemRoMedic™

WendyLett4Way

SystemRoMedic™

Instructions for use - English

WendyLett2Way and WendyLett4Way are two different draw sheets that can be used together with a WendyLett slide sheet. WendyLett2Way is used for lateral repositioning and WendyLett4Way is used for both lateral repositioning and repositioning higher up or lower down in bed. WendyLett4Way has an anti-friction material on the underside and has been developed mainly for heavy users. Available in several sizes.

Check Safety

Visual inspection

Check the condition and function of the assistive device regularly. Always inspect the product after laundering. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. If there are signs of wear, the product must be discarded.

Friction test


Low friction - Ensure that the surface is smooth and slides effortlessly. Compare with a new product.

Always read the manual

Always read the user manuals for all assistive devices used during a transfer. Keep the user manual where it is accessible to users of the product. Always make sure that you have the right version of the user manual. The most recent editions of user manuals are available for downloading from our website, www.directhealthcaregroup.com.


Placement

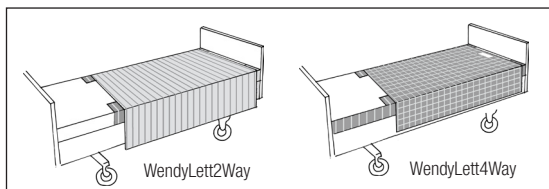Place WendyLett with the smooth surface in the middle of the bed and with WendyLett2Way on top. After turning, gliding movement can be inhibited by bedding WendyLett2Way in under the mattress on both sides of the bed.

A tag on WendyLett4Way indicates which side should be placed facing upwards in bed. After turning, gliding movement can be inhibited by bedding WendyLett4Way in under the mattress on both sides of the bed.

Material

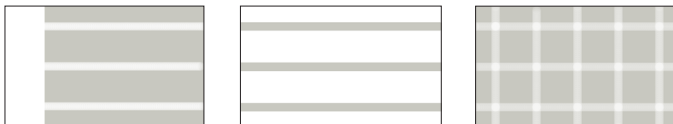The WendyLett slide sheet has a pattern of white stripes on a grey background. Stripes indicate direction of glide.

The draw sheet, WendyLett 2Way, has a white fabric with grey stripes. Stripes indicate the direction in which the sheet is drawn.

WendyLett 4Way has a grey-checked pattern and has a sliding surface underneath. Stripes indicate the direction in which the sheet is drawn.

Turning with WendyLett2Way


Perform the turn by grasping the draw sheet at the level of the hips and the shoulders. For best ergonomic technique, grasp a lot of fabric.

Work with straight wrists and shift weight by leaning backwards.

Turning with the aid of a draw sheet can be performed by two caregivers. This is a collaborative effort whereby one caregiver pulls the draw sheet (as above) and the other presses down somewhat near the shoulders and hips while at the same time pushing.

Turn the user. After turning, gliding is prevented by tucking the draw sheet under the mattress on both sides.


Turning with WendyLett4Way

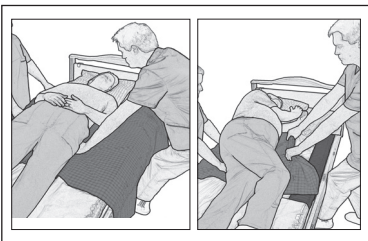The WendyLett4Way draw sheet, which has an anti-friction material on the underside, has been developed specially for heavy users. The procedure for turning is the same as for WendyLett2Way (above).

After turning, gliding is prevented by tucking the draw sheet under the mattress on both sides. To prevent chafing, smooth out the bedding under the user before tucking the draw sheet in under the mattress. The caregivers can do this by pulling diagonally on the draw sheet from either side. Raise the bed rails before leaving the user.


Higher up in bed with WendyLett4Way

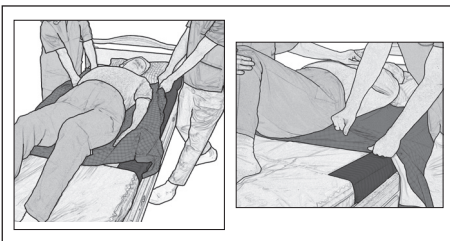To reposition the user higher up in bed, grasp the draw sheet on both sides near the shoulders, stand in a walking stance and take one step backwards.

To prevent the user from gliding down in bed, tuck in part of the draw sheet, so that it is folded double under the user's seat.




Sitting up in bed with WendyLett4Way

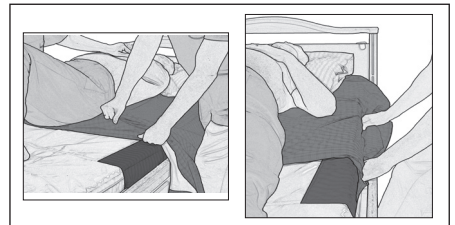With a WendyLett4Way in the bed, downward gliding can be prevented even when the user sits up in bed. Tuck in part of the draw sheet, so that it is folded double under the user's seat. Then, tuck the draw sheet in under the mattress, but only at the hip section. When the head end of the bed is subsequently raised, the draw sheet acts as a friction surface to prevent sliding.



Lift-assisted turning

When manual repositioning seems to place too much strain on the body, WendyLett4Way with handle straps (Art. No. 1657) can be used in combination with a user lift and a two-point slingbar.

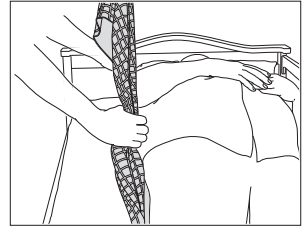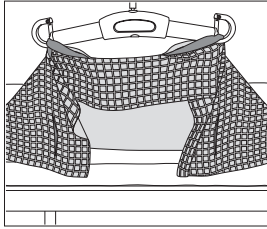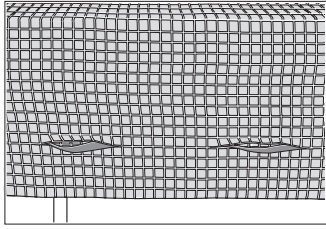Select two of the handle straps on one side and hook them to the slingbar. Prepare the user for turning.

Start turning carefully and stop to make any necessary adjustments; e.g., move the user closer to the middle of the bed to free more space for turning. Then, raise the slingbar until the user is lying on his or her side.

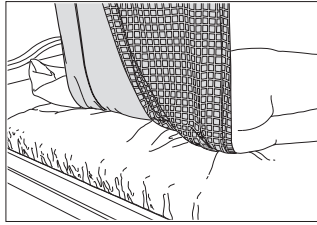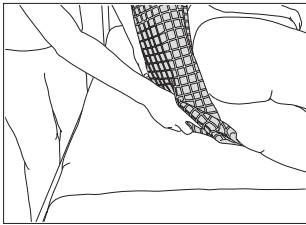With one hand under the hip, pull the lower edge of the WendyLett4way in, so that it is folded double.

The users back is now supported and the caregiver is easily able to attend to the users personal hygiene needs.

Note

Be aware that there is a risk of the user sliding out of bed. Raise the bed rails before leaving the user.

Material

Polyester, cotton, satin

Care of the product

Read the product label.

Do not use rinsing agent. To ensure maximum material life, hang dry or tumble dry on low heat.

Contact your local distributor if you have any questions about the product and its use.

See www.directhealthcaregroup.com for a complete list of distributors.

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Úplný zoznam distribútorov nájdete na www.directhealthcaregroup.com

Symbols

	This product complies with the requirements of the Medical Device Regulation 2017/745		Visual Inspection
	Medical Device		Read the manual
	Caution		Product Code
	Manufacturer information		Batch Code

Simple solutions for great results

SystemRoMedic® is the name of Direct Healthcare Group's unique easy transfer concept, the market's widest and most comprehensive range of clever, easy-to-use and safe transfer and lifting aids designed to make life easier, both for the user and for the caregiver. SystemRoMedic® is a complete solution that provides for the majority of patient transfer or manual handling requirements. From the simplest to the most complex scenarios, from the lightest to the heaviest. The concept encompasses assistive devices for four different categories of transfers:

- Transfer, assistive devices for manual transfers of users between two locations.
- Positioning, assistive devices for manual repositioning of users within the same location.
- Support, assistive devices for mobility support e.g., during sit-to-stand or gait training.
- Lifting, assistive devices for manual and mechanical lifting of users.

Improved work environment, improved quality of care and cost savings

The philosophy behind SystemRoMedic® is focused on the prevention and reduction of occupational injuries while allowing users to experience a greater sense of independence and dignity. Through a unique combination of training and a complete range of efficient transfer aids, SystemRoMedic® offers improvement of both work environment and quality of care and, at the same time, achieves significant cost savings.

Always make sure that you have the correct version of the manual

The most recent version of all manuals are available for downloading at/from our website;
www.directhealthcaregroup.com

For questions about the product and its use

Please contact your local Direct Healthcare Group and SystemRoMedic® representative. A complete list of all our partners with their contact details can be found on our website; www.directhealthcaregroup.com.

Moving Health Forward

Direct Healthcare Group Sverige AB
Torshamnsgatan 35,
SE-164 40 Kista, Sweden

Tel: +46 (0)8-557 62 200
info@directhealthcaregroup.com
www.directhealthcaregroup.com

Medical Device Class I. The product complies with the requirements of the Medical Devices Regulation 2017/745.