



# EASY HighBack

EASY™ in the high-back version (HB) is a development of the conventional concept using the back strap/neck support already available. This hygiene chair is fitted with a high, padded, water-repellent back complete with neck support (available in two versions). The high backrest provides soft yet stable support and its lower section is fitted with a pliant lumbar support for additional comfort. EASY™ HB is designed with telescoping pillars to create as much space and flexibility as possible for care personnel in the performance of their work. EASY™ is available with our new, soft Comfort seat with an opening at the rear for a more natural, hygienic working method.

**Integrated neck rest**

The integrated height-adjustable neck rest is available in two sizes, M and L.

**High, upholstered backrest**

The high backrest provides soft yet stable support. It is easy to clean and can be replaced as necessary.

**Tilt function**

EASY can be inclined 5 degrees forward and 30 degrees backward. This makes it easy to use when sitting down, getting up, and for personal hygiene. The tilt function also improves the working environment and reduces strain on employees.

**Pillars**

The telescoping aluminium pillars have an open design that allows water and dirt to simply run through.

**Swing-away, detachable footrests**

Flip-up, swing-away, detachable footrests with a heel restraint to prevent the patient's foot from sliding off in the tilt position are fitted as standard.

**Directionally locking castor**

The right front castor is lockable for easy manoeuvrability.



**Easy control**

Electrical height adjustment is via hardwired controller with a high ingress protection rating (IP65) for excellent safety.

**Lockable armrests**


The armrests can be flipped up or locked in their lower position. This is for safety when the chair is tilted.

**Pliant lumbar support**

The lower section of the backrest is fitted with a pliant lumbar support for additional comfort.

**Integrated seat with opening to the rear (Accessory)**

The seat's ergonomic design and flexible structure provide the best possible comfort for the user, while the opening to the rear provides the carer with a natural, more hygienic working method.

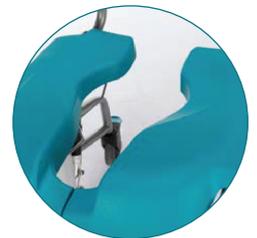

**Lockable castors**

All castors are lockable except the directionally lockable castor.

**Technical specifications**

Part No.	Product	Width, external	Max user weight	Seat height	Width between armrests	Seat depth
56-271	EASY HB	65cm	150kg	53-93cm	49cm	45cm
56-271/LOW	EASY HB Low	68cm	150kg	50-81cm	49cm	45cm

Examples of possible parts combinations are shown below.