

Thin, pliable and very easy to use

handicare

110 20

Perforated for the optimal format

An all-in-one solution for maximum hygiene and safety

## EasySlide on a roll – a new, unbelievably versatile sliding mat for short-term use

EasySlide disposable is a new, unbelievably versatile sliding mat that makes many common transfers, even heavy ones, easier for both the user and the caregiver. The practical, hygienic package contains a roll with up to 60 individual EasySlide mats for shortterm use. That means many easy transfers with maximum hygiene and safety.





# Easy transfers on a roll

EasySlide disposable is a thin, pliable and very easy-to-use sliding mat that can be used to facilitate turning and repositioning higher up in bed, for supine and seated transfers, and for repositioning in the seated position. EasySlide disposable is an unbelievably versatile assistive device that can also be used to lift the user's legs in and out of bed, to assist dressing and undressing, and to facilitate application of a lifting sling.

### A personal sliding mat for short-term use ensures maximum hygiene and safety

EasySlide disposable is ideal for use in hospitals, as well as in other situations and contexts in which extra high demands are placed on hygiene.

EasySlide disposable is made of a strong, environmentally friendly plastic fabric that has superior glide properties. It can be used as a single-use, disposable assistive device, but can also be used for short-term use with a single user. Instead of laundering, the sliding mat is discarded after use, when it becomes soiled or when the user no longer needs it. This prevents the spread of multi-resistant and other types of infectious bacteria while eliminating the inconvenience of washing with rubbing alcohol, growing piles of laundry, and the cost and problems associated with special laundering.

### Superior glide properties reduce friction and make transfers easier

EasySlide disposable is a tubular sliding mat, not a plastic bag. The outside has no discernible glide function; the sliding mat is easy to apply but the user needn't have to slide around uncontrollably. However, the inside of the sliding mat has lower friction, which gives superior gliding properties that allow the caregiver to perform safe and easy transfers, on virtually any surface, with a minimum of exertion and no unnecessary lifting. It is light and comfortable for the user, and easy on the caregiver's back and shoulders.

# Sliding mats on a roll – perforation gives the right format, every time

EasySlide disposable provides glide and easy transfers on a roll. Rolls come in two-packs, with each roll in a practical and hygienic plastic dispenser pack that can be hung on a wall or placed on tabletop or shelf. Each roll contains 60 slide mats, each 50 x 70 cm. Perforations at 50-cm intervals allow you to easily adapt the format of the mat according to the situation. Just unroll and tear off the desired length. The rest remains on the roll, clean and fresh, protected from wear and dirt, in the convenient dispenser pack.

#### Versatile and very easy to use

EasySlide disposable has many applications and is easy to use. The user guide is always at hand, since the instructions are printed directly on each mat. For turning and repositioning in bed, for supine and seated transfers, for repositioning in a seated position and for application of lifting slings

EasySlide

EasySlide

EasySlide

asySlide

EasySlide

18661

EasySlide

EasySlide

(COB BOAR

EasySlide

EasySlide



## Choose the format according to the current situation and needs

The smallest format, 50 x 70 cm, is ideal for heavy, relatively active users, and for positioning and repositioning a hip or shoulder in bed or on the x-ray or operating table. This format is also perfect for seated transfers, for example, from bed to wheelchair, shower chair or toilet, and for repositioning farther back in a chair.

With the 100 x 70 cm format you can reduce friction at both the shoulders and hips, thereby enabling quick and easy turning and supine transfers of heavy and somewhat passive users, for example, between bed and x-ray or operating table, shower trolley, stretcher, etc. Two 50-cm mats can be used to achieve the same effect.

With an even larger format, such as 150 x 70 cm or 200 x 70 cm, friction can be alleviated on all particularly sensitive points, including hips, shoulders and feet, when turning even heavier and completely passive users. The transfer can also be performed by sliding the entire mattress from a bed to, for example, an x-ray table.

EasySlide disposable is an assistive device that makes it easier for more active users to move on their own. Naturally, it's a matter of improving quality of life, but activating a user in bed also encourages circulation and micro-movements that prevent bedsores.

The sliding mat is made of a material that is not water soluble, making it the perfect solution for transfers in showering and hygiene situations.

# EasySlide, disposable



**Sliding mats on a roll.** Rolls come in packs of two, with up to 60 thin, pliable and easy-to-use sliding mats on each roll. Each roll comes in a practical and hygienic plastic dispenser pack that can be hung on the wall for safe, easily accessible storage.



EasySlide disposable is also available with a practical wall and tabletop dispenser in steel.

Art. No.	Product	Material	Description
2040 2045	EasySlide, disposable EasySlide, disposable		50 cm x 70 cm, 2 x 60 pcs 50 cm x 70 cm, 2 x 60 pcs, including practical wall and tabletop dispenser in steel

### CE 🛷 🖗 🖉 🖓 🔬

### Environmentally friendly material


The sliding mat is made of environmentally friendly polyethylene (PE) plastic, a durable material which, when incinerated, gives off only carbon dioxide and water. Used sliding mats can also be recycled as soft plastic. The packaging can be recycled as hard plastic.

### EasySlide disposable is an all-in-one solution:

- Always at hand
- Always new, clean and fresh
- Tough and durable, with superior glide properties
- Many different application areas
- Greater possibilities for activating the user
- Eliminates unnecessary and heavy lifting, which saves the caregiver's back and shoulders.



Perforated at 50-cm intervals; just unroll and tear off the desired length.

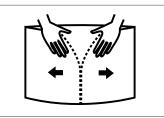

Tubular. Can also be used open; simply separate.

#### Other versions of EasySlide

EasySlide is available in several versions and sizes to suit different transfer situations and user needs. The classic versions of EasySlide are made of microfibre (polyester in 2010, 2011, 2015) with padding of polyester and ULF (Ultra Low Friction), our own, unique sliding material. Most sizes in these versions can also be fitted with a replaceable protective cover, which is made of polyester and can easily be cleaned. Disposable plastic covers provide for even better hygiene and less laundry.

Visit our website and read more about the other versions of EasySlide: www.handicare.com



SystemRoMedic<sup>™</sup> is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations. **Transfer**, between two locations. **Positioning**, within the same location. **Support**, for mobility. **Lifting**, both manual and mechanical.