

LECKEY^o

Stronger together.

We work with individuals, therapists and carers to design products with both a clinical and an emotional focus. Using the latest research and clinical understanding, we create practical solutions which are easily integrated into family life, because for us, life is about going, enjoying, participating and doing.



LECKEY^o

Leckey
Kilwee Business Park
Dunmurry
BT17 0HD
Northern Ireland
United Kingdom

T: 028 9060 2277
F: 028 9060 0795
E: info@leckey.com
W: www.leckey.com

LECKEY^o

Squiggles Saddle Seat

Squiggles

24 hour postural
care for babies,
kids & adults.
Sleeping, Sitting,
Standing, Walking,
Moving, Bathing,
Toileting.



The Squiggles Saddle seat has been designed to help develop active sitting. It promotes an abducted prone sitting position which helps achieve a neutral pelvic position providing improved trunk and head balance.


The combination of the contoured saddle seat, 4 point pelvic harness, adjustable footrest and back support provides appropriate and stable support enabling the child to focus on upper limb tasks.

The versatile Squiggles Saddle seat can be attached to either a Hi-low or Easy chassis and is ideal for one-to-one therapy work with a clinician or as a simple classroom seat for children with mild to moderate needs.

The backrest support has a unique design that rotates when the child fatigues so that you know when it's time for a break.

The colourful, tactile and fun design is ideal for young children and is available in 4 age appropriate colours.

Squiggles

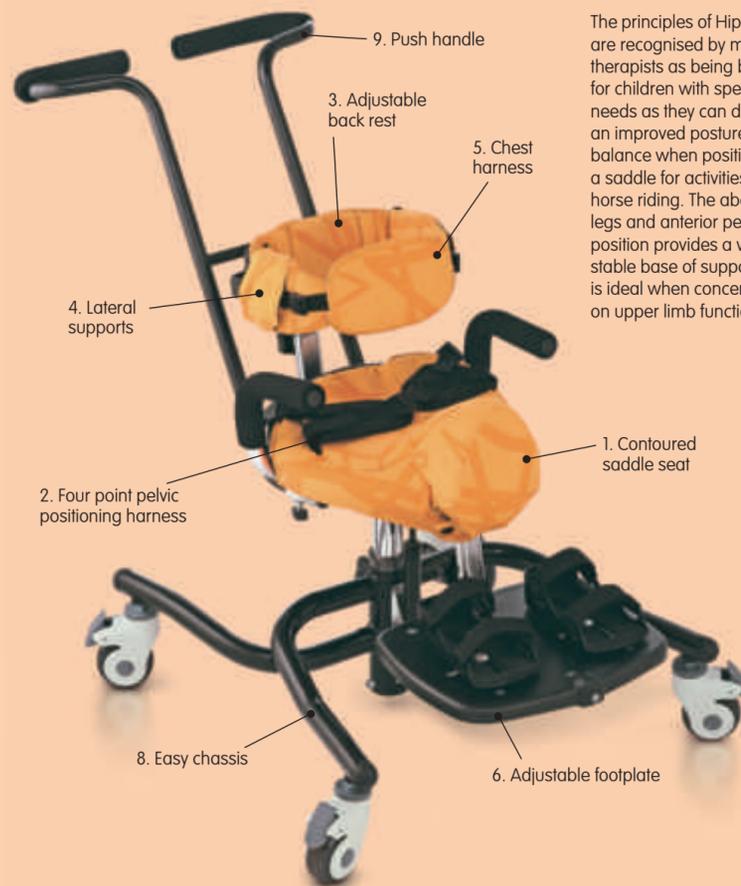

The principles of Hippotherapy are recognised by many therapists as being beneficial for children with special needs as they can display an improved posture and balance when positioned in a saddle for activities such as horse riding. The abducted legs and anterior pelvic position provides a wide, stable base of support which is ideal when concentrating on upper limb function.

1. Contoured saddle seat can be used to help develop active sitting or to provide an extra element of support when a child is concentrating in an active sitting position.
2. 4 point pelvic positioning harness holds the child comfortably in a secure but active position.
3. Back rest can be adjusted for the child's height.
4. Lateral supports help the child to maintain an upright position.
5. Chest harness also supports the child in the upright position.
6. Height adjustable footplate supports the child's feet and provides positive foot positioning.
7. The height and angle of the Hi-low chassis can be fully adjusted using a foot pedal for posture or activities.
8. The height and angle of the Easy chassis can be adjusted manually for posture or activities.
9. Push handle fits Easy chassis and Hi-low chassis.

Product sizes

Age (approx)	2 - 5	
Max user weight	22kg / 48.4lbs	
Stature	Min 930mm / 36.5 inches	Max 1105mm / 43.3 inches
Chest width (distance between laterals)	Min 160mm / 6.25 inches	Max 255mm / 10 inches
Seat depth (back support to back of knee)	Min 215mm / 8 inches	Max 290mm / 11 inches
Back support (Seat base to top of support)	Min 190mm / 7.5 inches	Max 265mm / 10.5 inches
Tilt in space (Easy chassis)	Prone 20°	Recline 20°
Tilt in space (Hi-low chassis)	Prone 10°	Recline 30°
Seat to footplate (Sole of foot to back of knee)	Min 190mm / 7 inches	Max 300mm / 12 inches
Footplate angle	Plantarflexion 10°	Dorsiflexion 10°
Armrest height (Seat base to elbow)	125mm / 5 inches	
Top of seat to floor (Easy chassis)	Min 280mm / 11 inches	Max 380mm / 15 inches
Top of seat to floor (Hi-low chassis)	Min 270mm / 10 inches	Max 650mm / 25 inches
Saddle seat unit (Min. dimensions for storage)	Weight 8kg / 17.6lbs Width 390mm / 15.4 inches	Length 450mm / 17.7 inches Height 370mm / 14.6 inches

The Squiggles Saddle Seat

The standard product includes:
Seat base; 4 point pelvic harness (small); Backrest; Chassis interface lock.

Seat Shell




122-600
Saddle Seat shell


Colour Options



122-751
Seat covers Green

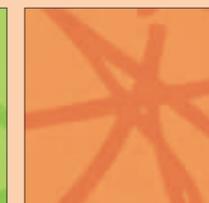

122-752
Seat covers Orange

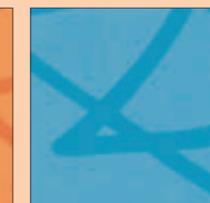

122-753
Seat covers Blue



122-754
Seat covers Pink

Covers include: support cushion and saddle seat cover

Accessories



117-781 Small sandals
117-782 Medium sandals



122-719 4 point pelvic harness



118-750 Support harness
+ 01 Green + 02 Orange
+ 03 Blue + 04 Pink



120-762 Rigid lateral supports
+ 01 Green + 02 Orange
+ 03 Blue + 04 Pink



122-664 Footplate assembly

Chassis Options

122-717
Easy chassis
manual
adjustment



120-800
Hi-low chassis
Foot pedal



122-710
Push bar handle
As shown with Easy Chassis
and Hi-low Chassis

Leckey Servicing

Our FREE service and support includes:
Free product training
Free product set up
Free product assessment
Free product re-assessment
Free repair within 2 year warranty

For service policies on all products outside warranty, please contact Leckey's Customer Service department. James Leckey Design Ltd as manufacturer with sole responsibility declares that all products conform to 93/42/EEC guidelines and EN12182 technical aids for disabled persons general requirements and test methods. Order forms and spare parts lists to extend the service life of the product and allow reissue are available on request or online at www.leckey.com.